

Dragon Medical One – Android PowerMic Mobile App Installation

For questions or assistance email speechrecognition@uhn.ca for all non-urgent requests.

PowerMic Mobile allows you to use your iPhone or Android smartphone as a secure wireless microphone for dictation with EPR. This application has been developed by Nuance for use with their desktop Dragon Medical One (Dragon) speech recognition application.

PowerMic Mobile for Android requires Android 4.0 or greater

If you do not have a QR scanner on your phone, open this installation guide on your smartphone. When you reach step 4 and 5, copy the URL below instead of using the QR code and skip to step 6 to complete the set-up.


http://config_/?NmsToken=MzlwMDMwRDktQUMwRi00RUZDLUEwQjMtMDA2MjEzMkJFNzAx&NmsBaseUrl=nms-ca.nuancehdp.com

Download and Install PowerMic Mobile


1. Open the Google Play Store on your Android and search for PowerMic Mobile. To download, tap **Install** located to the bottom of the application.



2. When ready, select **Accept** to agree to the end user license agreement.

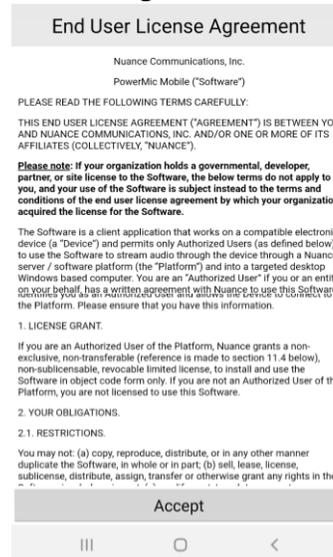



3. Select **Got it** from the Getting Started screen.

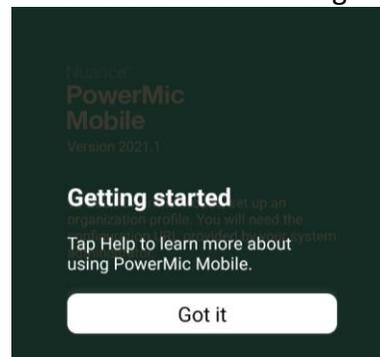


- Open the camera app on your phone and point it at this picture. If you have an older Android, press and hold the middle of the camera screen to activate the QR scanner.




- Tap the "Open in PowerMic Mobile" pop up if it appears or, tap Copy text to copy the configuration URL.

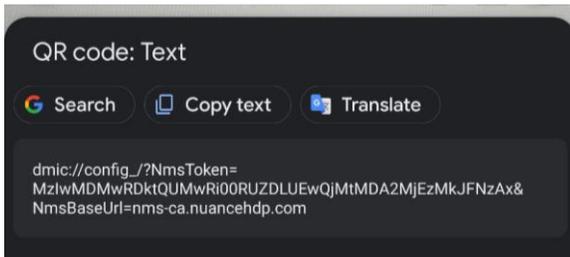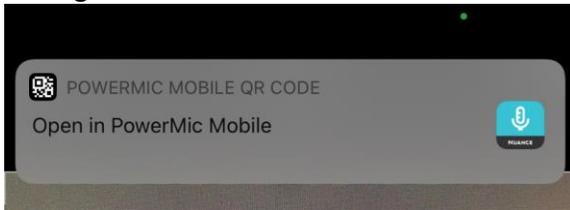


- If you copied the URL in the previous step, return to the PowerMic Mobile application and tap **Add Profile**

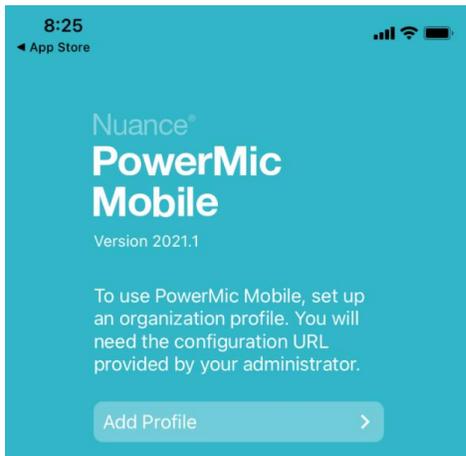


- Paste the URL into the **Profile URL** field.

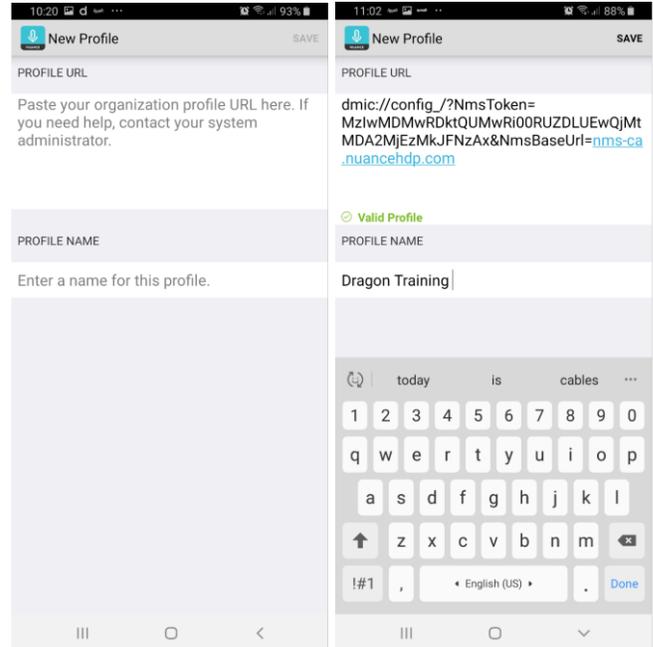


- Select **Save** located in the top right corner of your screen.

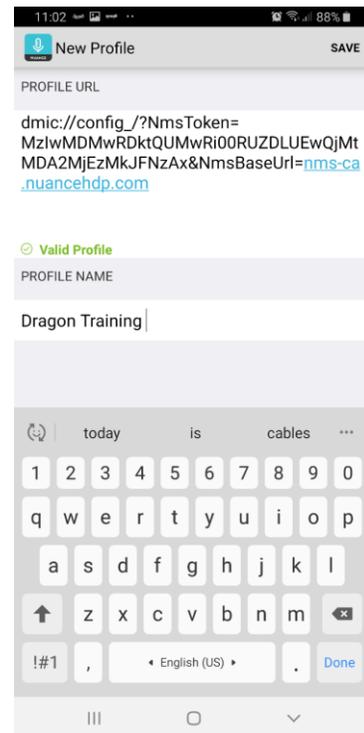


9. Tap Start to continue

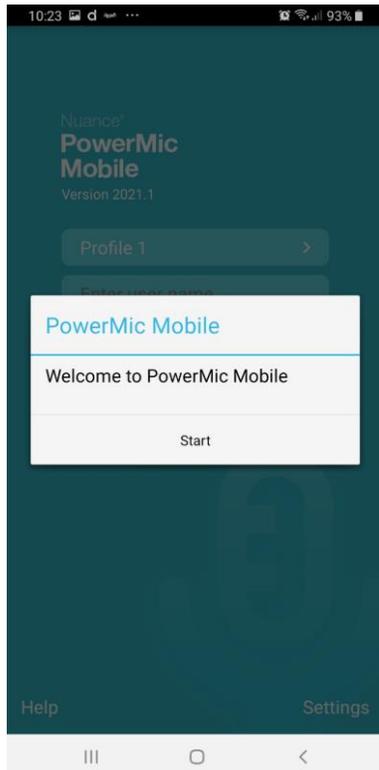


10. To ensure PowerMic Mobile Has Permission to Use Your Microphone, open the Settings and select Apps.

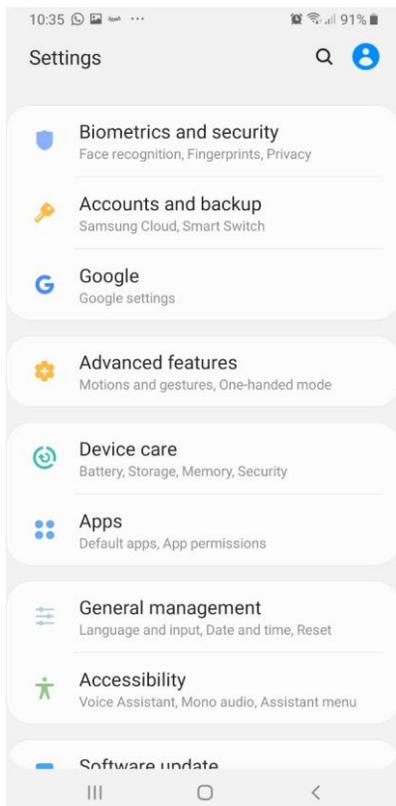

11. On the Microphone permissions page, ensure that PowerMic Mobile has access to your microphone. If it is greyed out, slide the switch to the on (blue) position.

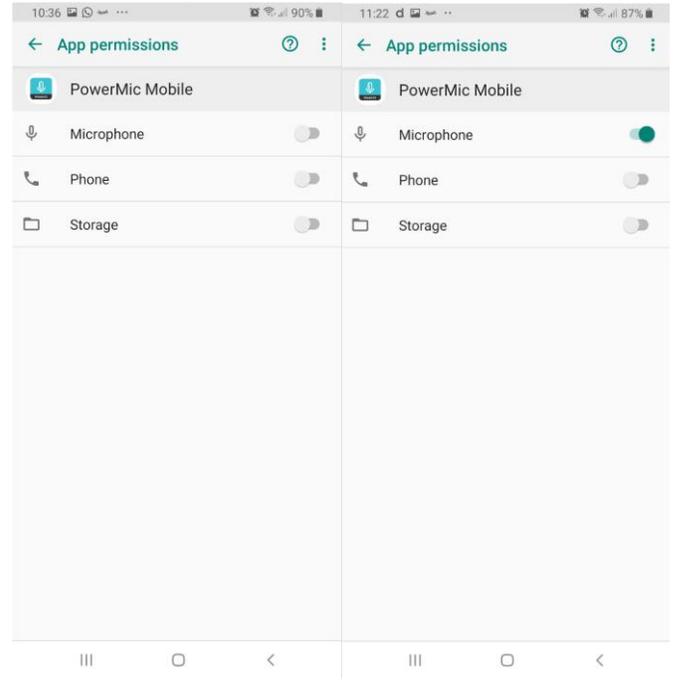


12. Log into PowerMic Mobile using your EPR ID as your username

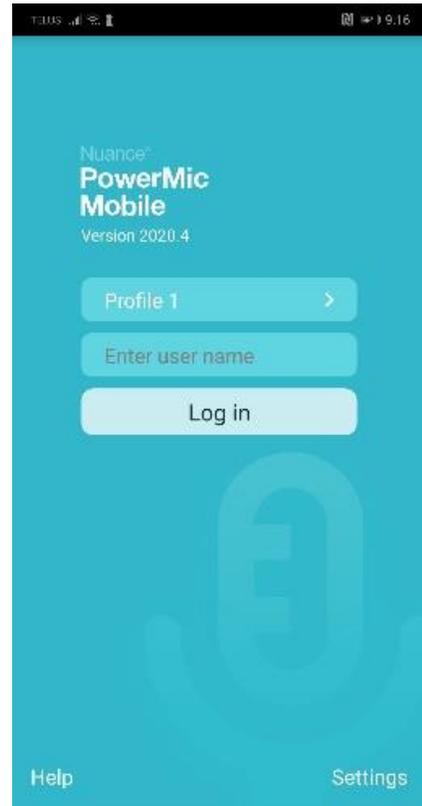

Pair PowerMic Mobile to Dragon Medical One

1. Launch the EPR with Speech application.

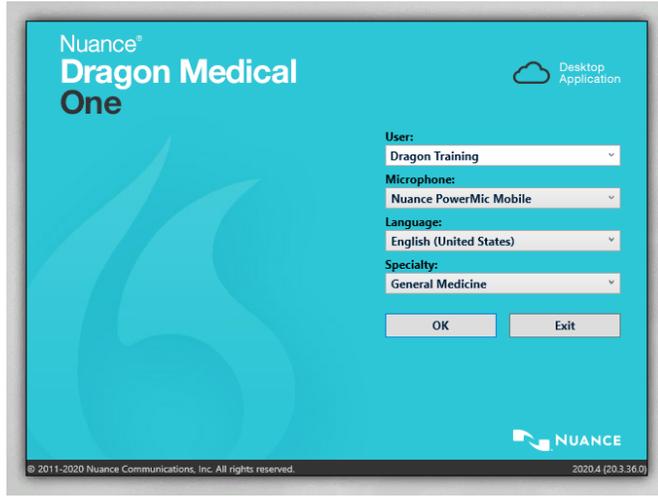


2. From the microphone field, select the drop-down arrow to select Nuance PowerMic Mobile and click **OK**.



3. With the PowerMic Mobile app open, pairing should be automatic. If the mobile phone icon on your Dragon bar has turned red, the pairing was successful.



4. If automatic pairing was not successful select the orange mobile phone icon on the Dragon bar to begin manually pairing your device.

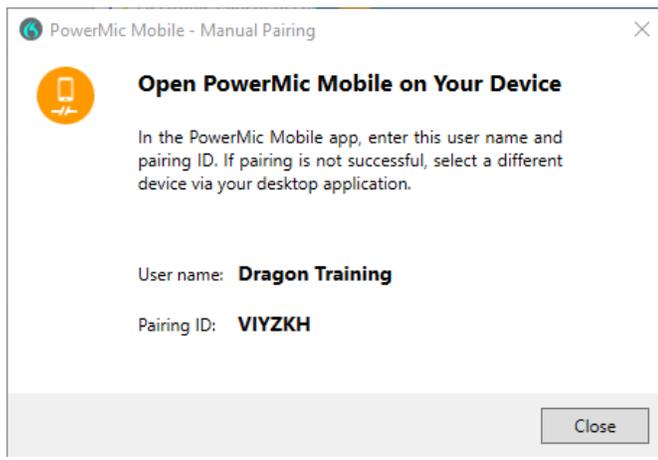

5. From the PowerMic Mobile application select Pair Manually and type the Pairing ID visible on your desktop.

